




COMPONENTS | DEVICES | POSSIBILITIES™

Urethane Marker Bands



SPECIFICATIONS

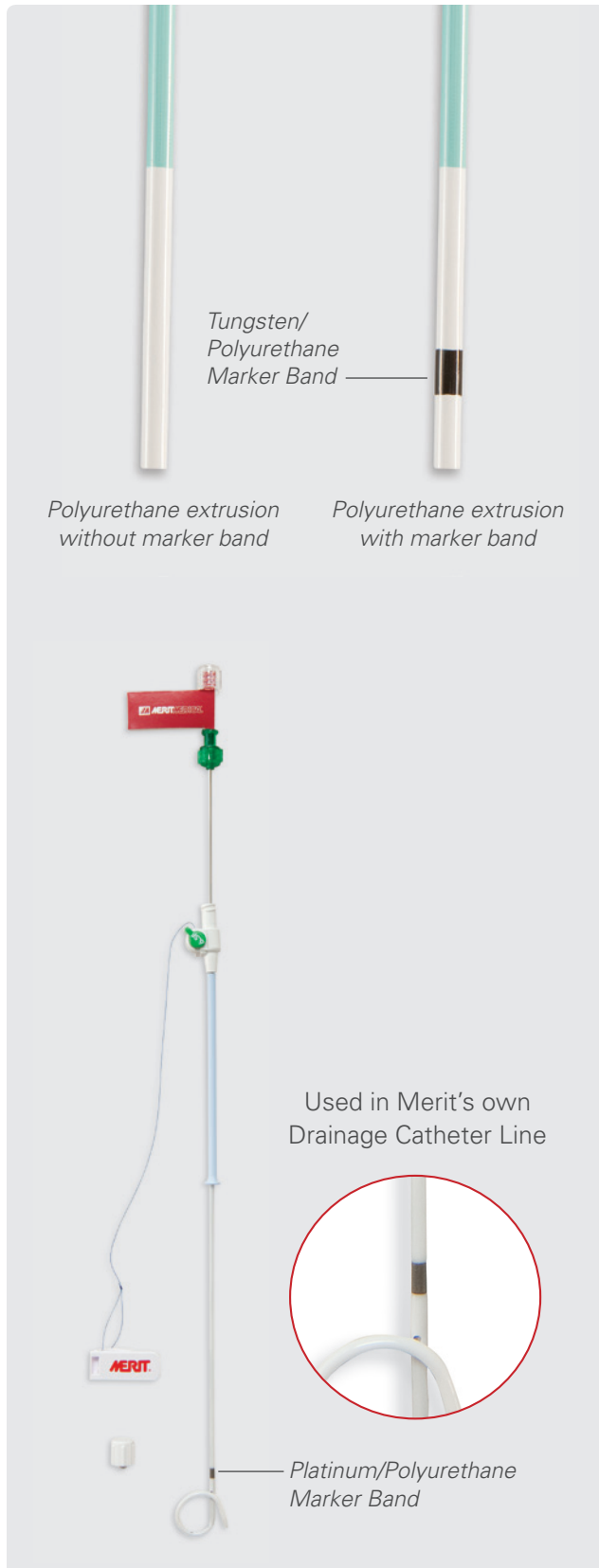

WHAT IT DOES

Marker bands help with guidance and placement of catheters during medical procedures. Merit's polyurethane marker bands, impregnated with tungsten or platinum, can be embedded into urethane family catheters. The marker bands are tough, flexible, and visible under fluoroscopy.

FEATURES & BENEFITS

Feature	Benefit
Embedded into catheter via melt-flow process; not crimped and swaged	Smooth transition; prevents kinking
Shaft OD's down to 0.020"	Can accommodate small OD sizes
Wall thickness from 0.002" to 0.006"+	Darken or lighten visibility of marker band under fluoroscopy

CUSTOMIZATION OPTIONS

Custom lengths and wall thicknesses available.

RELATED PRODUCTS

Interventional & Diagnostic Catheters



Drainage Catheters



HOW TO ORDER



Inside USA +1 800.637.4839
 Outside USA +1 801.208.4313
 Europe +31 43 358 82 24
 Or contact us at www.MeritOEM.com