


Inflation Devices

- A Worldwide Market Leader
- Ergonomic Handle Design
- Pressures up to 30ATM (440PSI)
- 5ml of safety space intended to minimize air bubbles introduced into the balloon
- Packaged stand alone or customized to your specifications
- Digital Technologies

May not be available for all OEM applications.

Single use. Sold packaged sterile. For non-sterile and other packaging options, please contact a Merit OEM representative.


SEE RELATED PRODUCTS

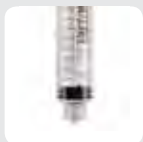

VacLok®
Syringes
Page 247



Stopcocks
Page 223

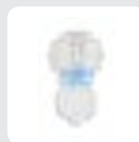

Hemostasis
Valves
Page 119

Monarch® Digital Inflation Device

Information at the touch of a button. A portable digital inflation device for the user who demands accuracy and convenience. Featuring selectable ATM/BAR or PSI readings and an easy-to-read backlit display. Available stand alone or call to customize to your specifications.

Display: Digital
Pressure: 30 ATM
Stopcock: 3-Way 500 PSI (34 Bar)
Tubing: 13" (33.02cm) Extension Tubing with Male Rotating Luer
Barrel
Volume: 20ml



- Timed Inflation Duration
- Timed Elapsed Since Last Inflation
- Displays pressures in ATM/BAR or PSI

Customization
Options Available

Monarch Inflation Device

IN2130

Blue Diamond™ Digital Inflation Device

Combining the accuracy of a digital device with the simplicity of an analog system. Featuring one button functionality for ease of use and a large dial face with oversize time and pressure digits to enhance visibility of critical information. Available stand alone or call to customize to your specifications.

Display: Digital
Pressure: 30 ATM
Stopcock: 3-Way 500 PSI (34 Bar)
Tubing: 13" (33.02cm) Extension Tubing with Male Rotating Luer
Barrel Volume: 20ml

NEW



- Timed Inflation Duration
- Timed Elapsed Since Last Inflation
- Displays pressures in ATM/BAR or PSI

Blue Diamond Inflation Device

IN7130

Basix™ COMPAK Inflation Device

Ergonomically designed to ensure accurate inflation of interventional devices. Features an angled analog gauge for better viewing of pressure readings and a luminescent dial for low-lit procedures. Available stand alone or call to customize to your specifications.

Display: Analog

Pressure: 30 ATM

Stopcock: 3-Way 500 PSI (34 Bar)

Tubing: 13" (33.02cm) Extension Tubing with Male Rotating Luer

Barrel

Volume: 20ml



Basix COMPAK Inflation Device

IN4130

Limited™ Inflation Device

No pressure display. Available stand alone or call to customize to your specifications.

Pressure: -1 to 25 ATM/BAR, 0 to 367 PSI

Tubing: 13" (33.02cm) Extension Tubing with Male Rotating Luer

Barrel
Volume: 20ml



Limited Inflation Device

K05-01299

Big60™ Inflation Device

Ergonomically designed to ensure accurate inflation and deflation of interventional devices. Features an angled analog gauge for better viewing of pressure readings and a luminescent dial for low-lit procedures. Designed to monitor inflation pressures over a range of -.68 to 12 ATM/BAR (-10 to 176 psi) with an accuracy of 3% ATM/BAR. Available stand alone or call to customize to your specifications.

Display: Analog

Pressure: 12 ATM/BAR

Stopcock: 3-Way

Tubing: 13" (33.02cm) Extension Tubing with Male Rotating Luer

Barrel
Volume: 60ml



BIG60 Inflation Device

ENDO-AN6012